



NEW PRODUCT GUIDE

Osteo-Arthritis Knee Brace #5539

The Osteo-Arthritis Knee brace is for mild to moderate unicompartmental knee osteoarthritis.

Features:

- Single side hinge provides 3 point pressure to shift joint forces from the degenerative area to the healthy compartment
- The hinge angle can be adjusted with the allen key provided
- Lightweight frame and anti-slip soft padding
- Quick release buckle for easy fastening and release

Universal sizing, specify left or right




Bort Asymmetric Knee #114700

Features:

- One lateral stay on each side provides lateral stability
- A silicone adhesive border provides a secure fit
- The DynamicFlex functional knitted fabric relax zone in the knee center relieves and protects flexor tendons and soft tissues
- Lengthwise elastic stretch zone in the lateral area of the silicone pad improves function and fit in all bending angles
- Moisture wicking fabric and antimicrobial finish wicks perspiration away, keeping the knee cool during activity and eliminating odor
- The top and bottom edges are made of extra soft snag proof material for extra comfort

Size	Calf Circ (d)	Thigh Circ (f)
S	26 - 29 cm	32 - 36 cm
M	29 - 33 cm	36 - 40 cm
L	33 - 37 cm	40 - 44 cm
XL	37 - 41 cm	44 - 48 cm
XXL	41 - 45 cm	52 - 56 cm

* Diagram applies to: 114420 / 114700 / 150140 -

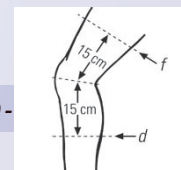

Bort StabiloGen Knee Support #114420

Features:

- Two lateral stays on each side provide lateral support
- The DynamicFlex functional knitted fabric relax zone at the knee center relieves and protects flexor tendons and soft tissues
- The lengthwise elastic stretch zone above the silicone pad ensures an ideal fit at all bending angles
- Moisture wicking fabric and antimicrobial finish wicks perspiration away, keeping the knee cool during activity and eliminating odor
- The top and bottom edges are made of extra soft snag proof material for extra comfort
- Compression class II

Size	Calf Circ (d)	Thigh Circ (f)
S	26 - 29 cm	32 - 36 cm
M	29 - 33 cm	36 - 40 cm
L	33 - 37 cm	40 - 44 cm
XL	37 - 41 cm	44 - 48 cm
XXL	41 - 45 cm	52 - 56 cm





Bort Stabilo Knee Support w/ Hinge #150140

Features:

- The extension limit is 0 degrees
- Circular fastening straps enable the hinges to provide extra lateral stabilization to the knee joint
- Pressure reduction zones in the extremities and popliteal fossa areas improve function and fit in all bending angles
- Moisture wicking fabric and antimicrobial finish wicks perspiration away, keeping the knee cool during activity and eliminating odor
- The top and bottom edges are made of extra soft snag proof material for extra comfort
- Compression class II

Size	Calf Circ (d)	Thigh Circ (f)
S	26 - 29 cm	32 - 36 cm
M	29 - 33 cm	36 - 40 cm
L	33 - 37 cm	40 - 44 cm
XL	37 - 41 cm	44 - 48 cm
XXL	41 - 45 cm	52 - 56 cm



Dynamic Knee Support

R5537

Features:

- 4 lateral stays for knee stability
- Gel patella pad for good patella support
- Cotton lining and breathable elastic material for comfort
- Contoured design for good fit and performance
- Velcro compatible nylon provides universal size

Sizing: Universal

Dynamic Wrist Lacer

R5593

Features:

- Lace closure system for optimal fit and easy application
- Breathable foam material for cool comfort
- Palm and dorsal stays for good support
- Reversible stays to fit either left or right
- Only 3 sizes to fit a wide range of wrist sizes

Size	Wrist Circumference
S	13 - 17 cm
M	17 - 21 cm
L	21 - 23 cm





Active Ankle Eclipse #77-1

Single upright multisport Ankle Brace. This ankle brace is engineered for maximum mobility. The innovative multipoint strapping system provides comfort and support to the ankle. Ideal for sports that require more mobility.

Size	Men's	Women's
S	5 - 8 1/2	6 1/2 - 9 1/2
M	9 - 12	10 - 13
L	12 1/2 - 16	13 1/2 - 16
XL	12 1/2 +	16 1/2 +

Active Ankle Eclipse 2 #77-2 \$39.90

Double upright multisport Ankle Brace. This ankle brace is engineered for maximum protection. The innovative strappings and padding system provides comfort and enhances the support from the double uprights. Ideal for sports that require mobility and protection against injury.

Size	Men's	Women's
S	5 - 7	6 - 9
M	7 1/2 - 9 1/2	9 1/2 - 12 1/2
L	10 - 12	13 - 16
XL	12 1/2 +	16 1/2 +



Dynamic LSO #5540

Features:

- Dual independent mechanical advantage pulley system for precise upper and lower compression, best fit, and maximum support
- Flexible pulley panels contour to the natural spinal lordosis
- Easy to use drawstring tabs for clients with poor hand strength
- Built-in abdominal panel for intra-abdominal pressure
- Removable pulley system for washing

Use product #5540B for black and #5540 for grey

Size	Waist Circumference
S	28 - 34"
M	34 - 40"
L	40 - 46"
XL	46 - 52"

Dynamic LSO High Profile #5542

Features:

- Precontoured posterior panel provides optimal support and pain relief
- Posterior panel is made from heat moldable ABS to custom mold if required for a custom fit
- Posterior panel is covered with a breathable liner for comfort
- Modular posterior panel is removable for a lower profile fit over time

Sizing: See Dynamic LSO sizing guide





Edge LSO

#5561

LSO anterior height: 9"
LSO LP anterior height: 7.5"
LSO & LSO Low profile
posterior height: 14.5"

Features

- Rigid spinal system to provide maximum stability
- Intracavitary pressure to unload intervertebral discs
- Back panel opening for wound inspection
- Comfortable moisture wicking padding
- Low profile version also available (product #5560)

Size	Waist Circumference	Hip Circumference
S	30 - 36"	35 - 41"
M	34 - 40"	39 - 45"
L	38 - 44"	43 - 49"
XL	42 - 48"	47 - 53"



Edge TLSO

#5561T

Features

- Rigid spinal system to provide maximum stability
- Intracavitary pressure to unload intervertebral discs
- Back panel opening for wound inspection
- Comfortable moisture wicking padding
- Low profile version also available (product #5560T)

Size	Waist Circumference	Hip Circumference
S	30 - 36"	35 - 41"
M	34 - 40"	39 - 45"
L	38 - 44"	43 - 49"
XL	42 - 48"	47 - 53"

Hallux Valgus Night Splint

#5580

Features

- High quality splint for night time correction of Hallux Valgus
- Fabric lined padding for durable, comfortable support
- Easy to apply velcro strap

Universal sizing, specify left or right



Cool Exchange Compression Bandage #80695 / #80696

80695 (Regular) - Recommended for ankles, feet, wrist, elbows and fingers - 3" x 47" - \$16.00

80696 (Large) - Recommended for knees, thighs, groin, hamstrings, back, shoulder or cut to fit - 3" x 118" - \$32.00

CoolXChange is an innovation gel bandage which advances R.I.C.E. (Rest, Ice, Compress, Elevation) methodology by combining the cooling and compression stages, providing a convenient and effective 2-in-1 solution for pain relief. Heat energy is transferred out via evaporation so no refrigeration is required. Compression and cooling minimizes swelling. Can be worn continuously and comfortably for extended periods without ice burns. To reuse simply put it back in the package with 3 teaspoons of water. CoolXChange will rehydrate for up to 20 more uses.